DE/EN

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**INSTRUCTIONS
FOR USE**

AssistX ABL one

CSS
MICRO SYSTEMS

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WEEE DE94635635

Table of Contents

1. Introduction	18
2. Intended use	18
3. Explanation of symbols	19
4. Safety instructions	20
5. Cleaning and disinfection	23
6. Scope of delivery	24
7. Device Buttons, Ports, Display	25
8. Use of the device	26
9. Settings	26
10. Maintenance / Review	27
11. Support	27
12. Technical specifications	28
13. Conformity / Standards	28
14. Environmental protection & disposal	29

1. Introduction

Thank you for choosing a high quality CSS product. Please read these instructions for use carefully and take care of these instructions when using the device. Please keep the instructions for use at hand when you use the device.

2. Intended use

The **AssistX ABL one** is a reassurance lamp for the call system interface AssistX Call. It is used to supplement call triggers that do not have their own reassurance lamp (e.g. proximity sensors, special buttons).

The **AssistX ABL one** is connected to the **AssistX Call** and placed within the caller's field of vision. If a call is triggered via the **AssistX Call** and acknowledged by the nurse call system, the red LED light of the **AssistX ABL one** lights up and a sound comes up (can be switched off).

3. Explanation of symbols

Safety instructions!

Important instructions!

4. Safety instructions

Handling

Please handle the device carefully. The device contains sensitive electronic components. It is not allowed to do any modification. The device is not life-sustaining or life-supporting. The device consists of several components. Small items may be detached from the device, which may result in suffocation or other danger to children. Small children or people with intellectual disabilities should not reach or use the device without supervision.

Do not use the device in damp, wet, very cold or very hot environments (see temperature information in technical data).

Never replace the built-in batteries by yourself! The battery can be damaged, which can lead to overheating and can lead to injuries. The integrated rechargeable lithium polymer battery may only be replaced by authorized personal.

Mounting

The device must be mounted in a way that it cannot fall down. Only recommended mounting components from the manufacturer are to be used. The user bears the sole risk of mounting.

The manufacturer is not liable for personal injury or property damage caused by the device or its components falling down.

Repair

Under no circumstances you should open the device and repair it by yourself. Disassembling the device can permanently damage the device and cause injury. Contact your reseller or the manufacturer if the device is damaged or does not work properly. If the device or its accessories are mechanically damaged, they may no longer be used.

Interaction with medical devices

The device and its components can emit electromagnetic radiation and have integrated magnets. The device must therefore not be placed on a person or brought closer than 15 cm to people using other active medical devices. The device must also not be placed directly on active medical devices. Please pay special attention to people with implanted pacemakers. If necessary, please contact the device manufacturer of the other active medical device.

Attention! Magnets

The mounting plate contains a strong neodym magnet. Magnets can affect cardiac pacemakers and implanted defibrillators.

- The pacemaker could be switched into testmode and cause discomfort
- A defibrillator could stop working
- Keep enough distance to magnets if you use such devices
- Advice other users of such device against getting too close to magnets
- Keep magnets away from all devices and objects, which can be damaged by strong magnetic fields

Recommended safety distance is 6 inch (15 cm).

Magnets cause a far-reaching magnetic field. They could damage TV sets, laptops, PC harddrives, creditcards, ATM cards, data storage, mechanical watches, hearing aids and speakers.

Vehicle keys, USB sticks, cameras, cell phones and smart phones contain non-magnetic storage media. Therefore, static magnetic fields near those devices cannot delete data. There are countless covers with magnetic closures or magnetic holders on the market, which also speaks against the notion that magnets can damage these devices. However, it cannot be ruled out that very strong magnetic fields might magnetise and maybe damage mechanical parts or built-in speakers. When in doubt, keep these devices away from strong magnets.

5. Cleaning and disinfection

Before cleaning the device, switch it off and disconnect all cables. The device should be cleaned if required with a disinfectant and cleanser intended for use in hospitals. Cleaning should only be done with a wipe. Wipe all surfaces and let them dry or follow the wiping instructions.


Do not use sprays that may allow moisture to enter the device. Under no circumstances may any component of the device be subjected to steam or wet sterilization.

6. Scope of delivery

- **AssistX ABL one**
- 2 m 4pol. connection cable
- Instructions for Use

7. Device Buttons, Ports, Display

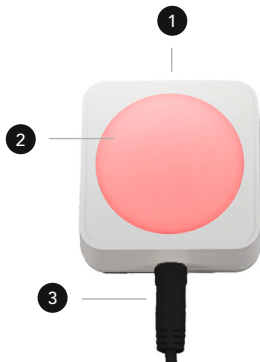- 1 On and Off Button (Sound)
- 2 Red confirmation LED
- 3 Input (**AssistX Port**)

8. Use of the device

The **ABL one** is connected to one of the **AssistX Port** connections (S1 or S2) on the **AssistX Call** using the supplied cable. The device is then immediately ready for operation.

The **ABL one** lights up and an audio signal sounds (can be switched off) when the call system has confirmed that the call has been triggered.

9. Settings

The confirmation sound can be switched on and off on the **ABL one**. (default: switched on) By pressing the ① button for more than 5 sec, the sound is switched off (the ② LED lights up briefly once, beeps once). To switch on, press the ① button for more than 5 sec: the sound is switched on (the ② LED lights up 2x and beeps 2x).

10. Maintenance / Review

The following checks must be carried out every three months:

- Function
- Connector for damage and dirt
- Damage of the housing and other components

11. Support

If you need any support or for technical questions about the product, please contact the reseller where you purchased this product. If necessary, e. g. the reseller is not available or cannot help, send us an email to support@csslabs.de. Please be sure to include a description of the problem and the serial number of the device.

12. Technical specifications

Dimensions: 50 x 50 x 26 mm

Weight: 60 g

IP Rating: IP 54

Power supply: 12VDC via **AssistX Port**

Connection: **AssistX Port**, 3.5 mm jack, 4pin cable

Temperature ranges

Use: 0 °C to 45 °C

Transportation and storage: -10 °C to 50 °C


13. Conformity / Standards

 This product is CE compliant.

AssistX ABL one meets the requirements of DIN VDE 0834 in conjunction with **AssistX Call** and **AssistX Mobil**.

14. Environmental protection & disposal

All electrical and electronic devices must not be disposed with general household waste, but must be disposed at designated disposal points. This product is subject to the European Directive 2002/96/EC and must therefore be properly disposed at the WEEE collection points.

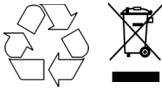WEEE DE94635635

THE PASSIONATE ASSISTIVE TECHNOLOGY MANUFACTURER

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© 2022 CSS MicroSystems GmbH